

Sediment

(Marine Sediment, PACS-1)

Sample Category: Environmental & Regulatory

Sample Requirements

Sample Weight: 0.5g

Reagents: 10 mL Conc. HNO₃

Vessels: 35 mL glass, quartz, and Teflon™ lined

Microwave Heating Program

Control Style: Ramp To Temperature


Pressure Mode: Organic

Recommended Equipment



Discover SP-D

Stage	Temp (°C)	Ramp Time	Hold Time	Pressure (psi)	Power (W)	Stirring
1	200	04:00	02:00	400	300	Med

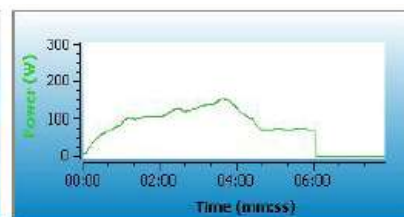

Notes:

1. If using a vial other than the recommended choice, adjust sample size and pressure limit to values appropriate for the vessel chosen.
2. This method may not provide a total digest of all sediment samples.

General Precautions:

This procedure is a reference for sample digestion using a CEM System and may need to be modified or changed to obtain the required results on your sample