

Vitamin

(Ultra Mega Gold)

Sample Category: Food / Plant / Animal Tissue

Sample Requirements

Sample Weight: 0.25g

Reagents: 3mL Conc. HNO₃, 9mL Conc. HCl

Vessels: 35mL glass, quartz, and Teflon™ lined

Microwave Heating Program

Control Style: Ramp To Temperature




Pressure Mode: Organic

Recommended Equipment



Discover SP-D

Stage	Temp (°C)	Ramp Time	Hold Time	Pressure (psi)	Power (W)	Stirring
1	200	5:00	3:00	400	300	Med

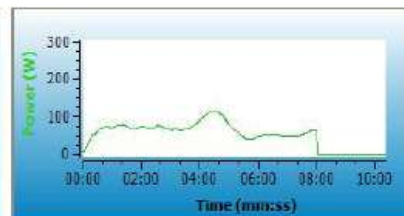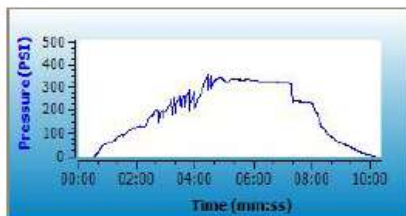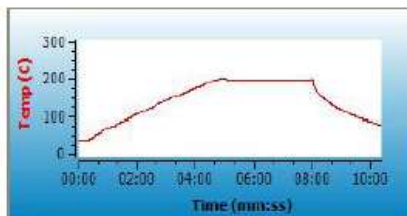

Notes:

1. If using a vial other than the recommended choice, adjust sample size and pressure limit to values appropriate for the vessel chosen.

General Precautions:

This procedure is a reference for sample digestion using a CEM System and may need to be modified or changed to obtain the required results on your sample